

How to install MQLink


Symbion Lavery Pathology uses the latest in message delivery technology for its Electronic Data Interchange (EDI) purposes. MQLink is a point to point message delivery system utilising IBM WebSphere and TCP/IP Internet protocols. MQLink can deliver messages in many data formats and is user configurable for many of its options. Each individual site that receives data is allocated a unique site code and all data is encrypted before transfer and decrypted once delivered to the site computer.

System Requirements:

- Windows NT, XP, 2000 or 2003 Operating System
 - * Please contact EDI Support if running Windows Terminal Server or ISA Firewall.
- Minimum 20 Mbyte disk space
- Pentium III processor or better
- Existing Internet connection
 - *MQLink will not initiate an Internet connection
- Existing clinical software
 - *Only required for processing of results, not result delivery itself.

Installing MQLink:

1. Ensure you are logged in as Administrator, or, at least have Administrator rights.
2. If not already done, Download the mqlink.exe file from the Symbion website and save to your hard drive (or run from supplied EDI disk).
3. Double click the mqlink.exe file.
4. Answer Yes to "Install Symbion Health Results Delivery System".
5. The installation will unzip and commence.
6. Enter your *site code* and *check digit*.
 - * If you do not already have this information phone (02) 9005-7819 for details.
7. In the Results Folder, enter the required destination of the transferred file.
 - * Mapped drives and/or UNC pathnames are acceptable.
8. If your site will be receiving Diagnostic Imaging results as well as Pathology, enter the required destination in the DI Results Folder.

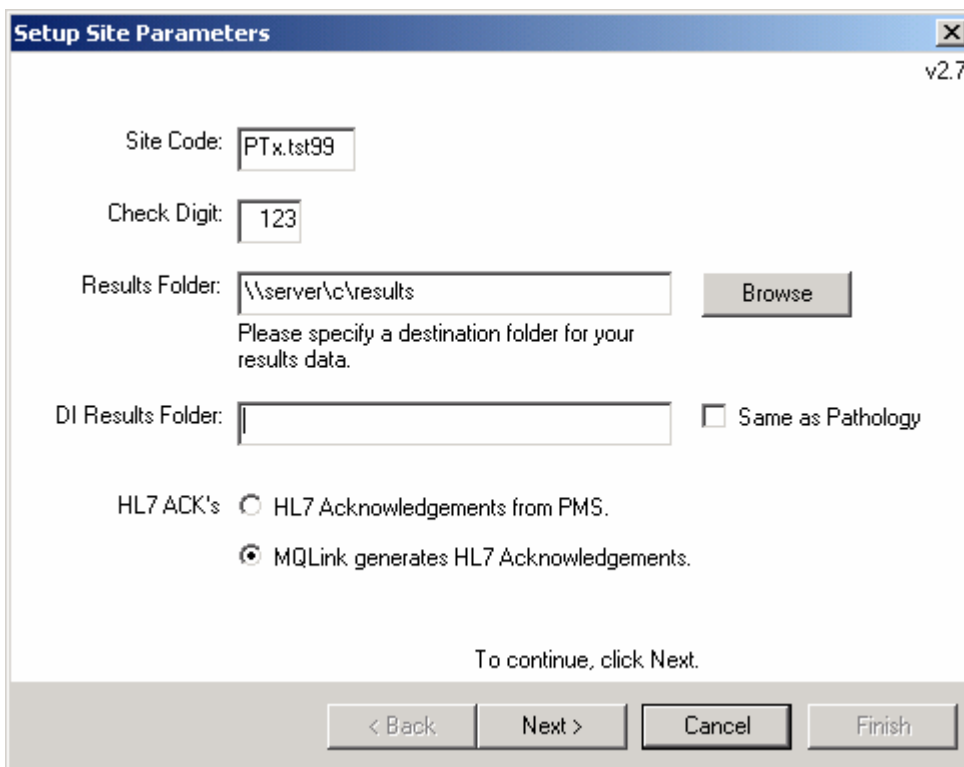


(**'Same as Pathology'** tickbox will duplicate the Results Folder information).

9. Leave the HL7 ACKS's setting as is unless directed by Symbion EDI staff.
10. Click Next.
11. When the install has completed, click Finish.

****NOTE**** MQLink is a single installation system. **DO NOT INSTALL ON MULTIPLE COMPUTERS!**

Starting MQLink:



1. Click on Start > Programs > Startup > MQLinkStartup.
The MQLink Pathology Results Transfer System Icon will appear in the system tray. *(Bottom right hand corner of your screen next to the system clock)*

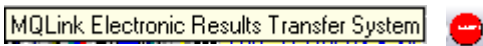

2. Right click on this icon then left click on '**Collect Pathology Results**' to commence the download.

***Note** The first run of MQLink will generate a site-specific encryption key and upload it to our servers. You will not receive any results during this initial run).*

The task will appear on your task bar. Click this button to watch the progress.

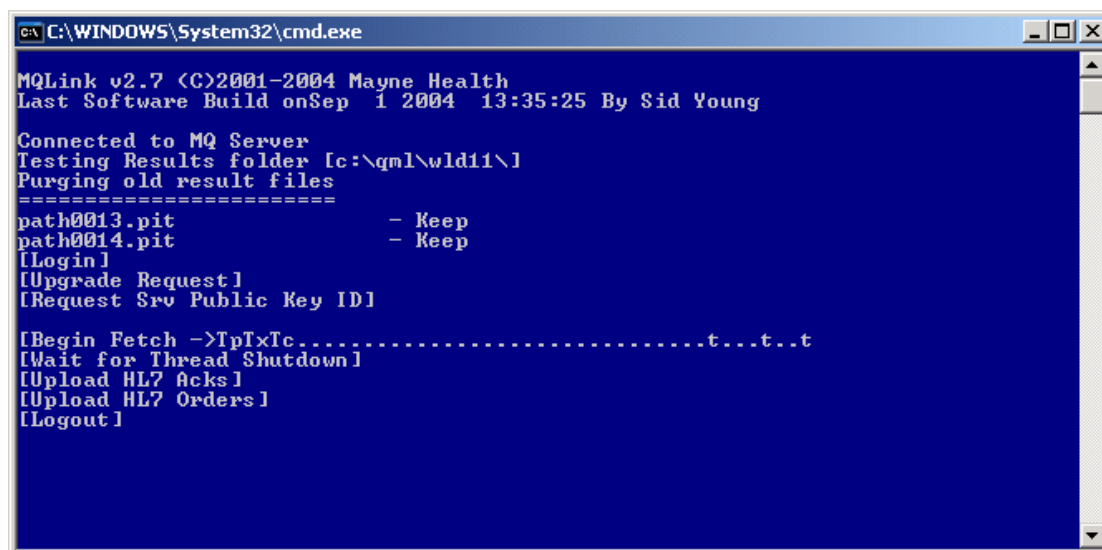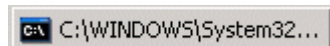

When complete, the window will close.

3. To schedule the downloads to occur automatically, right click on the MQLink icon and left click on '**Task Scheduler**'. Enter the details for the Start Time, End Time and required Interval times of the downloads. Click on '**Set Event**' to save the information.
4. The Download installation is complete and will now run automatically.

Contact Symbion Lavery Pathology EDI staff on (02) 9005-7819 to arrange for release of data or for technical assistance.